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**Sterilization Autoclave Class B 18L / 23L
with a build-in printer**

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Contents


1. Introduction	29
2. Safety	29
3. Packing list and installation	32
4. Description and specifications	33
5. Control Panel	36
6. Operation.....	37
7. Logs Reading software.....	45
8. Maintenance	48
9. Transportation and storage.....	52

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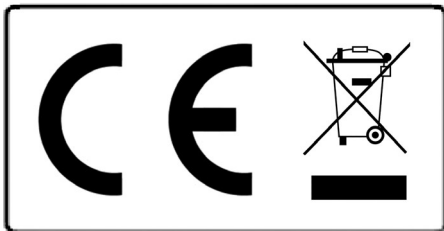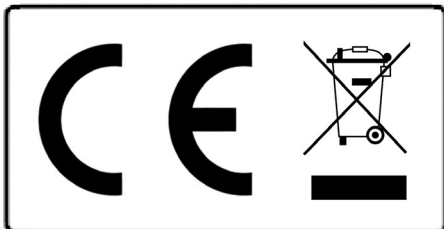
Notice to customers:

The product is made up of various components and various materials that must be recycled or, failing that, deposited in the corresponding debris removal sites when the product's life has been completed or when otherwise it is necessary to dispose of it. To do this, the end user who acquires the product must know the current regulations of each municipality and / or locality based on the waste electrical and electronic equipment. The user who acquires this product must be aware of and responsible for the potential effects of the components on the environment and human health as a result of the presence of hazardous substances. Never place the product in a conventional container of citizen scope if a previous dismantling and knowledge of the components that incorporates. If you do not know the procedure to follow, consult with the city council for more information.

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1. Introduction

1.1 Important

- . This manual contains necessary and sufficient information to operate the system safely, like optimal use, a safe and reliable operation, regular and correct servicing requirements.
- . Read and understand all instructions in this manual before attempting to use the product.
- . Keep this manual with the equipment at all times. Periodically review the procedures for operation and safety precautions.

1.2 Usage indications

This equipment can be used for dental clinic, laboratory, surgical room, emergency room, ophthalmology, gynecology and steam, cosmetic hospital, etc ... and also used by doctors and professionals.

1.3 Contraindication

There is not any contraindication of this equipment.

2. Safety

2.1 Safety symbols

Device icon description

	ATTENTION: refer to this operating manual. It alerts the user to refer to the operation manual or other instructions when complete information cannot be provided on the label.
	ATTENTION: High temperatures - Pay attention to high temperatures inside the chamber and outside the sterilizer when the exhaust system is working.
	Protective Earth - It indicates the protective earth terminal (ground).
	CAUTION: Hazardous voltage (arrowhead beam) used to indicate electric shock hazards.

Label icon description

	Symbol for "SERIAL NUMBER"
	Symbol for "CATALOGUE NUMBER"
	Symbol for "DATE OF MANUFACTURE"
	Symbol for "MANUFACTURER"
	Symbol for "AUTHORISED REPRESENTATIVE IN THE EUROPEAN COMMUNITY"
	Symbol for "CAUTION"

Manual Prompt

Note	Indicates that concerning information which easier or helpful in operation
Caution	Indicates that a potential hazard may exist which through inappropriate conditions or actions will or can cause: <ul style="list-style-type: none"> • Minor injury • Property damage. • Damage to machine
Warning	Indicates that a specific hazard is known to exist which through inappropriate conditions or actions may cause: <ul style="list-style-type: none"> • Severe personal injury • Substantial property damage • Substantial damage to machine

NOTE: Indicates precautions or recommendations that should be used in the operation.

2.2 General Safety Recommendations

- The user is responsible for operation and servicing the sterilizer in accordance with the instructions listed in this manual.
- The sterilizer could not be used for liquid.
- The sterilizer has not been designed to operate in the presence of gas or explosive vapors.
- The trays and the load will still be hot at the end of each cycle. Use the tray holder to remove each tray from the chamber.
- Do not open the door of the chamber during the sterilization programs.
- Do not put your hands or face on the cover of the water tank when the sterilizer is running.
- Do not remove the instruction plate or any label from the sterilizer.
- Do not pour water or any other liquid over the sterilizer.
- Do not fill the caustic liquid into the water tank.
- Do not fill the caustic matter in the chamber.
- Use only high quality distilled water.
- Unplug the mains lead before inspecting or servicing the machine.
- Only an approved technician using original spare parts can carry out repairs and maintenance.
- In case of transport, drain both water tanks completely, allow the sterilization chamber to cool down and preferably use the original packaging.
- The articles under sterilization should be removed by special tools when the temperature over 40°C.
- Picking-up the sterilizing trays should use the special tools provided.
- During the transportation, the sterilizer should be carried by two people in case of turning over.
- Notice! This product can't be put on the place that is not easy to cut off power supply.
- Prohibit covering the cover of water tank during usage.

2.3 Safety parts and functions

Temperature Protection:

- Temperature Protector (Steam Generator) Cut off current when the steam generator temperature is too high.
- Temperature Protector (Heating Ring) Cut off current when the heating ring temperature is too high.

Electronic Protection:

- Double Fuse It will be broken circuit if the connected power is too high or unstable.
- Electronic Filter It will filter the electromagnetic interference during working.

Mechanical Protection:

- Jiggle Switch To ensure the door is closed completely for avoiding the risk of safety.
- Tray Tong Avoid scald when removing articles from the chamber.

Control Part:

- Temperature Sensor (Inside Chamber) To measure temperature inside the chamber.
- Temperature Sensor (Heating Ring) To measure temperature of the heating ring.
- Temperature Sensor (Steam Generator) To measure temperature of the steam generator.
- Pressure Sensor To measure pressure of the chamber.
- PCB Control Control system for controlling all the process of sterilization.

CAUTION!! MANUFACTURER IS NOT HELD RESPONSIBLE FOR ANY ARBITRARY DISASSEMBLY, AMENDMENT WITH THE UNIT, BY UNAUTHORIZED PERSON OR UNPROFESSIONAL TECHNICIAN.

2.4 Operation risks

Please take attention to avoid the following operation risk during operation:

Scald Risk

- Every time open the door after finish the cycle, please keep an appropriate distance, because the chamber still have rudimental steam with high temperature, avoid to be scalded.
- Every time open the door after finish cycle, please not touch the main door and chamber, because of high temperature, and avoid to be scalded.

Pollution Risk

Please clean the chamber after every time use to avoid rudimental contamination left inside chamber.

2.5 Protective device

Plastic or cloth gloves are protective devices useful during loading and removing articles, to avoid scald.

3. Packing list and installation

3.1 Packing list and accessories

The standard equipment consists of the following components. Open the package and check whether everything is in good conditions. Keep the original package for some days, just in case.

Components	Pcs
- Wire tray	3
- Tray support	1
- Handle for handing trays	1
- Water drainage tube	1
- Water filling hose	1
- Electrical mains connection cable	1
- USB Pendrive	1
- Door gasket (spare part)	1
- T-Bubble Levelling device	1
- Water filters	1+1
- Fuse Holder	2
- User manual M80416	1
- EU Declaration of Conformity	1


3.2 Location

Place the machine on a flat, leveled surface, suitable for the equipment weight. Use the supplied leveller to level the autoclave.

Note: The autoclave should be slightly raised from the front side to ensure the water in the boiler not to accumulate back in the drain area. Otherwise, the «drying» phase will be incorrect and may cause the error message E3.

You will need a plug, no longer than 1m, in accordance with the equipment features. There's no need of a water tap nor a drain.

To work comfortably, leave about 50cm around the equipment. To ensure air circulation, leave a free space about 20cm behind the equipment. Be careful not to block any ventilation grid.

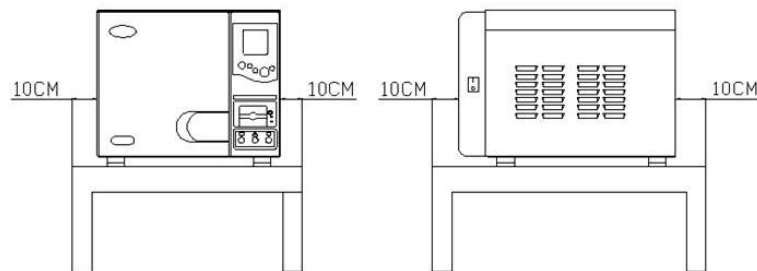The sterilizer should be set in a place which at least has 10cm distance with each side (20cm of the top) and should be set in a place with good ventilation.

The temperature of environment: 5-40°C.
 The humidity of environment: ≤85%
 Atmosphere pressure: 860Hpa-1060Hpa
 An earth connection is essential.

CAUTION: DO NOT PUT ANY STUFF WHICH EASILY MELTS NEAR THE STERILIZER.

- The sterilizer should be set on level table or place; the front-end should be a little bit higher than back-end.
- The sterilizer cooling and vent area should not be jammed or blocked.
- Do not put any stuff on the top of the sterilizer.
- Do not put any stuff in front of the door, to avoid accident when open the door.
- Do not put any corrosive stuff near the sterilizer to avoid accident or risk.

3.3 Connection to the mains

Locate an outlet (type 16A) near the autoclave. Use the network connection cable supplied or a similar one (section 1.5mm²)

The plug must be grounded.

Do not handle power cords with wet hands.

Check that the plug voltage matches with the indicated on the features label.

Connect the cable to the back side of the autoclave.

CAUTION:

- DO NOT BEND THE POWER CABLE TO AVOID DAMAGE
- DO NOT PUT ANY HEAVY STUFF ON THE POWER CORD
- DO NOT USE OTHER POWER CORD TO AVOID DAMAGE OF STERILIZER
- DO NOT TRY TO ADD THE POWER CORD TO AVOID ACCIDENT OR RISKS

3.4 Filling the water tank

Use deionized water for the autoclave. Fill the tank through the top hole. (3,5L for the 18L autoclave and 4L for the 23L one).

Another filling method is by means of the water inlet 3. For this, connect the provided hose immerse the other end into a container with fresh water. Press key 2 for 3 seconds to start the filling process. The equipment stops the process when it is completely filled.

3.5 Draining

The autoclave has two separate tanks, one for clean water and the other for used water. Both tanks are separately emptied, through the valves on the right front side. Connect a hose to open the valve.

4. Description and specifications

Autoclaves Class B are widely used in medical, dental, veterinary and generally to sterilize materials that may come in contact with blood, as for instance, tattoo centres, piercing salons, hairdressers or beauty centres, etc.

Autoclaves Class B assure proper sterilization of different types of normal processes: solid, solid packed, wrapped porous and hollow. To ensure sterilization, Autester ST DRY PV autoclaves achieve a complete air removal and a better steam penetration thanks to an efficient system of fractioned vacuum.

AUTESTER ST DRY PV Clase B print are a series of autoclaves which have the following features:

- Sterilization chamber and cover in stainless steel.
- External epoxy coated case.
- External steam generator to the sterilization chamber.
- Operation digital control.
- Specific programs for testing and autoclave diagnosis.
- Large backlight LCD display.
- Air inlet filter to the camera.
- Water tank.
- RS-232 output for printer and USB connection for logs downloading.
- Build-in printer.

To guarantee sterilization of the load, the appropriate program indicated in the instructions must be used.

Front View


Name

- 1. Water Filling Port
- 2. Drainage Port
- 3. Water Filling Port
- 4. Drainage Port
- 5. Internal Printer
- 6. Power Supply Switch
- 7. Displaying Window
- 8. Door Handle

Description

Manually filling water
 Connected to used water tank
 Automatically filling water
 Connected to main water tank
 Internal printer to print sterilization record after finished cycle.
 Standard green power switch
 It will show temperature, pressure or other sterilization information on the window, easy for you to operate during the working.
 Door handle with safe door locker to open and close the door.

Rear View

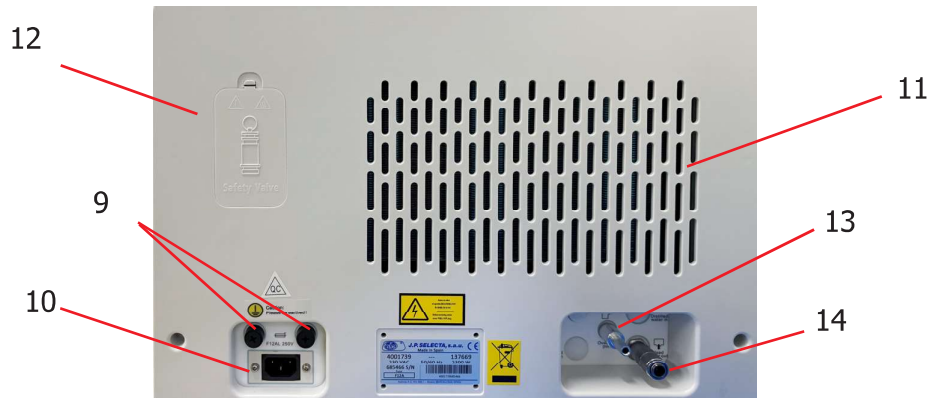Name

- 9. Fuses
- 10. Power socket
- 11. Vent area
- 12. Safety Valve
- 13. overflow connection
- 14. Drain Port

Description

Protect device from overcharge or electrical problem.
 Connected with power source
 Output heat from this vent area by condenser
 Protect device from overpressure
 To discharge water from used water tank if not manually empty by means of port (2)
 To empty used water tank, it's connected to port (2)

Open view

Name

- 15. Air Filter
- 16. Door Gasket
- 17. Tray Frame

Description

Filter the air and ensure air into chamber is clean
 For sealing the door
 To load the articles

4.1 Autoclaves type «B»

Autester ST DRY PV autoclaves allow sterilization of:


- Solid loads resistant to temperature: instrumental, metallic.
- Solid loads sensitive to temperature: plastics
- Solid packed load.
- Porous load: textiles
- Wrapped porous load.
- Catheters, tubes and cavities devices.
- Liquids.

4.2 Technical specifications

Model	4001738	4001739
Supply voltage:	220V, 50Hz	220V, 50Hz
Power:	1800VA	2300VA
Fuse:	F12AL	F12AL
Working temperature:	5-40°C	5-40°C
Noise:	< 50dBA	< 50dBA
Maximum weight per tray:	1 Kg	1 Kg
Chamber material:	Stainless steel	Stainless steel
Maximum working pressure:	2.5barg (250kPa)	2.5barg (250kPa)
Minimum working pressure:	-0.9barg (-90kPa)	-0.9barg (-90kPa)
Maximum temperature:	145°C	145°C
Chamber diameter measurements:	245x360mm	245x470mm
Water volume per cycle:	0.3L	0.35L
Water tank volume:	3.5L	4L

4.3 Sterilization cycle

- 1-2 pre-heating
- 2-3 pre-vacuum
- 3-4 sterilizing
- 4-5 air-discharging
- 5-6 drying
- 6-7 stabilizing
- 1-7 entire duration

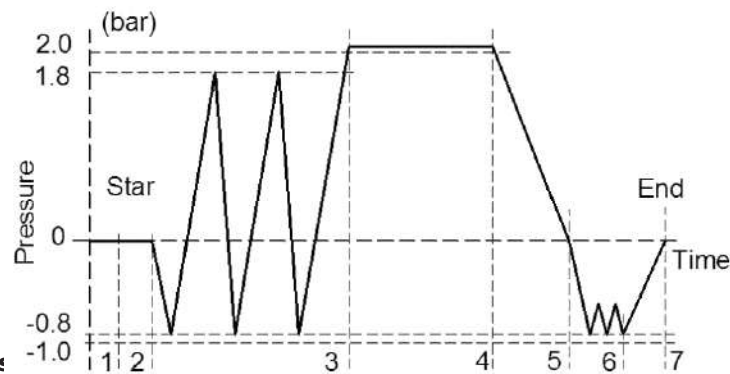Sterilization cycles type:

Type B: The sterilization of all wrapped or non-wrapped, solid, hollow load products type A and porous products as represented by the test loads in this standard.

Type S: The sterilization of products as specified by the manufacturer of the sterilizer including non-wrapped solid products and at least one of the following: porous products, small porous items, hollow load products type A, hollow load products B, single wrapped products, multiple-layer wrapped products.

NOTE 1 The description identifies ranges of products and test loads.

NOTE 2 Non wrapped sterilized instruments are intended either for immediate use or for non sterile storage, transport and application (e.g. to prevent cross infection).

4.4 Preset programs

Programs (Cycles)

Program	1	2	3	4	5
Program type	Unpacked	Packed	Prion	Rubber	Textile
Sterilization temp	134°C	134°C	134°C	121°C	121°C
Sterilization pressure	204kPa	204kPa	204kPa	104kPa	104kPa
Sterilization	4 min	5 min	18 min	20 min	20 min
Drying	9 min	9 min	9 min	9 min	18 min
Vacuum pulses	1 vacuum	3 vacuums	3 vacuums	1 vacuums	3 vacuums

Program	6	7	8	9
Program type	Liquid	User	B&D Test / Helix	Vacuum Test
Sterilization temp	121°C	134°C	134°C	
Sterilization pressure	104kPa	204kPa	204kPa	
Sterilization	30 min	4 min	3,5 min	
Drying	1 min	3 min	3 min	
Vacuum pulses	N/A	3 vacuums	3 vacuums	5 min - 5 min - 10 min

*Also apply to Helix Test.

5. Control Panel

- 1 Display
- 2 Increase value key
- 3 Decrease value key
- 4 Program key
- 5 Start/Stop key

5.1 Controlling Panel

5.1.1 Integrated display

- A. Display the chamber pressure, unit is: bar
- B. Display the chamber temperature, unit is: °C
- C. Display the cycle situations. Please see the circumstance which is attached (Sterilizer operation cycle display chart).
- D. Display the alarm error code when the sterilization gives the alarm during the working cycle, then you can check and find the malfunction components according to the error code.

5.1.2 Increase value key (2)

UP choosing has the functions that are used for setting down the cycle program and filling in water pump, Press this button you can choose the following programs, and if press button for 3 seconds, the machine would enter into filling in water automatically program, and the LCD would display the working cycle state.

5.1.3 Decrease value key (3)

DOWN choosing has the cycle program, press this button you can also choose the following programs.

5.1.4 Program key (4)

Press this button for 3 seconds, it will enter into the cycle program choice menu, and then please press (1) or (2) buttons to choose the cycle program as you need.

5.1.5 Start/Stop key (5)

This button is used for confirmation on the start or end the cycle and removing the malfunction of program. After you choosing a certain program, please press this button to confirm it. Press the machine again will start a working cycle. Keep pressing it for 5 seconds during a program any time, the cycle will be terminated.

CAUTION: PLEASE CONTACT YOUR DISTRIBUTOR OR AUTHORIZED PERSON WHEN ERROR CODE IS DISPLAYED.

6. Operation


6.1 Starting

Switching On

Please connect the power cord before you start the operation. The power switch to control the unit turn on or off. It is located back right corner of the sterilizer. When you turn the power on, you may see followed window:

Door is open before turn on power switch	Door is closed before turn on power switch
Close the door and then do next step	Open the door first, then to close the door follow the information on window, and do next step

Distilled Water Adding

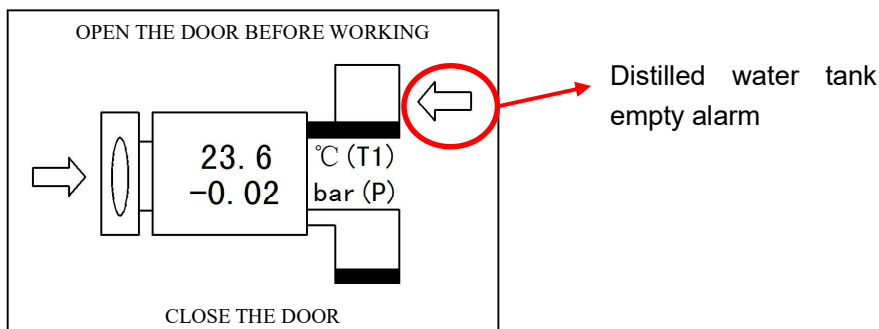After you connecting to the power and select the corresponding program, if the screen shows an WATER FILLING ARROW as above, which means the water in the water storage tank is of its lowest level and you have to add water in.

If the ARROW appears on the screen, even you press "start/stop" button, it will not work either, so you need to add distilled water in until the lowest required level can be reached. You can fill in water automatically or by manual on the top of the machine.

Automatically: Connect the water tube with the machine left lower corner inlet connector as right picture arrow marked, please see the picture, and press UP button for 3 seconds; by manual: open up the top cover, to fill in water according to the pic 3, stop filling in water again when you hear beep sound.


If the water shortage alarmed during the sterilizing, you should not worry about it, the water will not influence this process. The water that had been left in the tank is adequate for the operating cycle, but you need to add water in time, so as to ensure the next operation.

CAUTION: USE DISTILLED WATER ONLY TO EXTEND THE LIFE TIME OF THE MACHINE. DO NOT TILT THE STERILIZER WHEN THE RESERVOIR IS FULL OF WATER.

Alarming if Used Water Reservoir is full

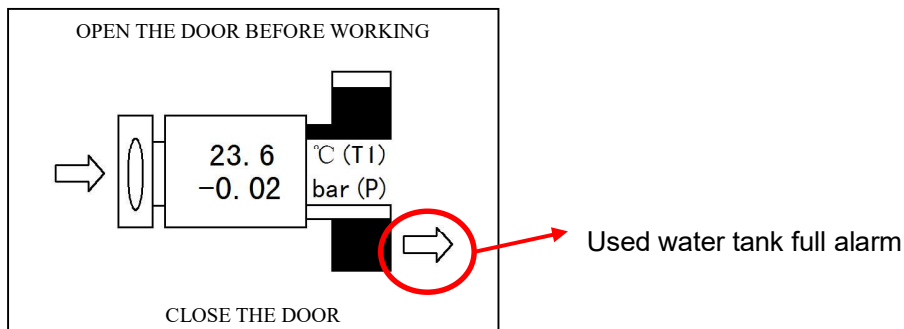The LCD screen shows WATER DRAINING ARROW during the sterilizing, which means that the used water reservoir needs to be drained out.

Connect the water tube to the machine left lower corner inlet connector as right picture arrow marked, please see picture.


Used water outlet

Generally, the maximum temperature of drained water should be under 70°C. If it is higher, you need to check whether the fan operates normally, or contact the local distributor immediately, we will offer our best service to you in no time.

6.2 Program selection

- Turn on main switch.
- Press PROG key (4) to access the options PROGRAM and CONFIGURATION.
- Press the UP key (2) to access the preset programs.
- Select the program with keys (2) and (3).
- Confirm by pressing PROG key (4).


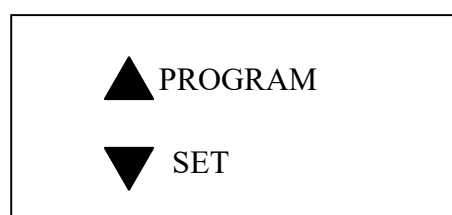
6.3 Sterilization

- Load the trays with the material to be sterilized.
- Press Start/Stop key (5) to start the selected program.

If this operation is made with the door closed, the display will ask to open the door and close it again before pressing the Start/ Stop key.

- The display will indicate the process status every moment.
- At the end, END will appear. Wait for the 0 bar condition.

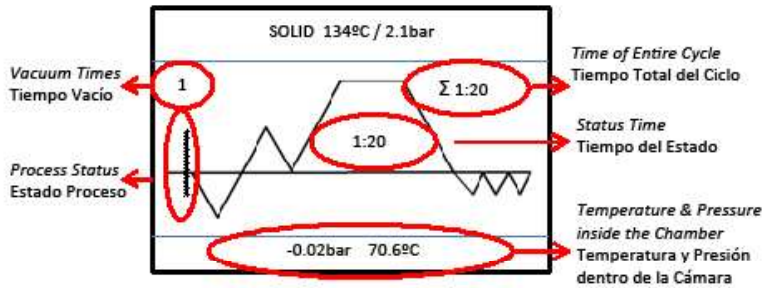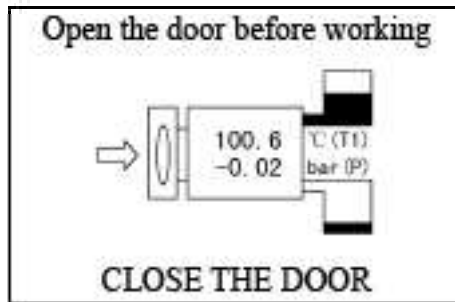Process graph that appears on the display


End of cycle

Some examples of PACKED program screens:

Program example with 3 pre-vacuum: PACKED

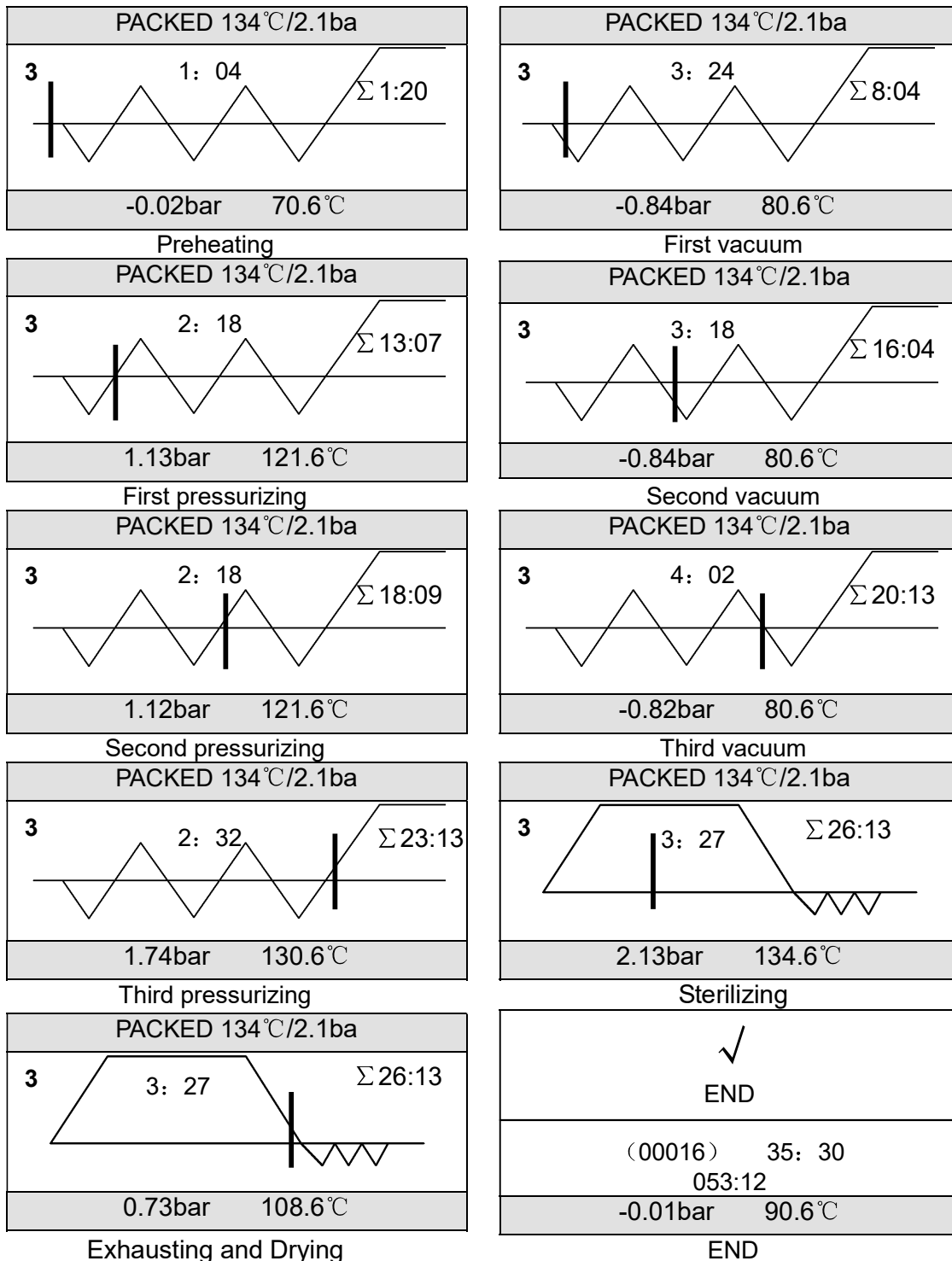Force end of program

In the cycle, if the program is interrupted by the error or pressing Start/Stop button (5) more than 3 seconds it will enter "Abnormal Exiting program" as followed window:


WARNING	
Abnormally exit	
0.73bar	108.6°C

In this state the air release solenoid valve will be opened and exhaust air. You have to cancel this information by pressing Start/Stop button (5) for more than 3 seconds to leave and to go back to normal state window.

WARNING: Do not try to open the door when the pressure value is higher than 0.1 bar.

6.4 Configuration

- Turn on the main switch.
- Press the UP key (2) to access the configurable parameters (Time, Printer, Language and Water Quality).
- Select the option to modify with keys (2) and (3).
- Press Start/Stop key (5) to access the submenu or choose the desired option.

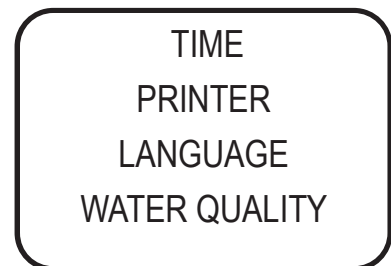6.4.1. Set time

D	/	M	/	Y
12		05		2012
H	:	M	:	S
11		02		20

Press UP or DOWN button to increase or decrease values, then press PROG button to change items, and the press START/STOP to save and exit.

6.4.2. Set printer

PRINTER
ALL / ON / OFF

Press UP or DOWN button to turn ON or OFF the printer, choose ALL to print the date which record the value in each 30 secs.
 OFF - neither the printer nor the USB works
 ON - the printer and USB works
 ALL - is factory only

6.4.3 Set language

LANGUAGE	
1. English	русский
2. Español	6. Italian
3. Português	7. French
4. 中文	8. Polski

Press UP or DOWN button to choose language, and then press START/STOP to save and exit.

6.4.4 Program 11 USER settings

11	User-defined	11	User-defined
	134°C 4min 03min 3time		▶ 134°C 4min 03min 3time
	↑ ↓ SELECT		↑ ↓ SELECT

Sterilization temperature	121°C	134°C
Sterilization time	15 to 50 minutes	4 to 40 minutes
Drying time	1 to 30 minutes	1 to 30 minutes
Vacuum number	0, 1 or 3	0, 1 or 3

Program No. 11 "User-defined" is used to achieve a perfect sterilization by adjusting to a special setting. The process of the "User-defined" program is as follows:

- Press the PROG key (4) to enter program 11, then press PROG key (4) again to select the value you wish to modify: sterilization temperature, sterilization time, drying time and number of vacuums.
- Press the UP (2) or DOWN (3) keys to change the value of each parameter.
- Press the START key (5) to confirm and START (5) again to start the cycle.

CAUTION when removing the material from the autoclave.

Risk of burns due to contact with hot surfaces.

The material inside the autoclave could be at a temperature higher than 100°C.

Use the handle provided for trays handling.

6.5 Loading articles


Instruments should be put on the instrument trays with some gap between each articles so that the steam can be ventilated freely. Please use the attached handle tool to load instrument trays into chamber to avoid scalding.

Arrangement on Trays before Sterilization

- Read the following instructions for proper usage and maintenance of articles and material.
- Make sure that the articles of different materials are separated and placed on different trays.
- In case of carbon steel articles, place a towel or paper-wrap between the tray and the articles in order to avoid a direct contact.
- All the articles must be sterilized in an open position.
- Make sure that the articles remain apart during the sterilization cycle.
- Do not overload the trays.

Drawings explain:

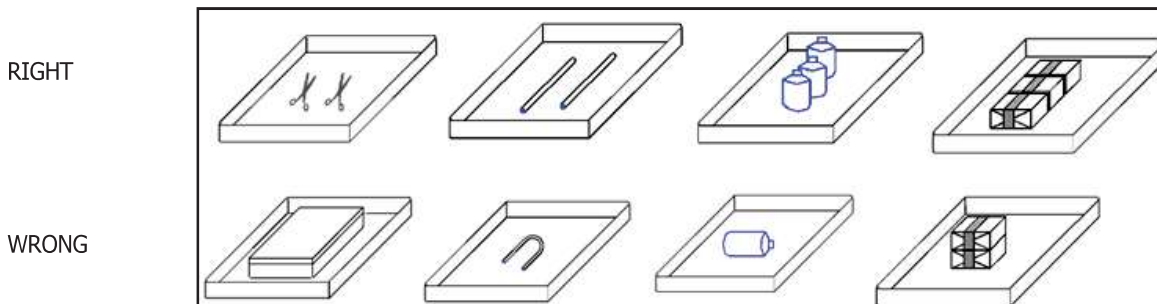CAUTION: RECOMMENDED TO CLEAN THE ARTICLES BEFORE LOADING.

ATTENTION: IT IS HELPFUL TO TURN ON POWER FOR WARMING UP FROM 5 TO 10 MINS BEFORE RUNNING PROGRAM, IF THE ENVIROMENT TEMPERATURE IS UNDER 10°C.

6.6 Closing the door

Close the door after loading the articles. LCD will not display alarm when latch the door handle fully. When close the door, however, if the chamber is warm and steam still be left in it, you may feel a strong resistance when you close the door. You just have to push harder and latch the handle completely. You may also leave the door open to release the steam and kindly close the door again. Or you can push the door in while you turn the door handle. Anyway, latch the door handle completely. If you still unsure the door has been closed properly. You can also adjust the door

If the door not close well, the LCD screen will show "close the door" and remind you to close the door completely.


CAUTION: THE DOOR SHOULD BE SHUT TIGHTLY BEFORE RUNNING PROGRAM TO AVOID THE DANGER.

ATTENTION: THE LCD WILL DISPLAY AN ALARM IF THE DOOR IS NOT CLOSED FULLY, THE STERILIZER WILL NOT RUN UNLESS THE DOOR HAS BEEN CLOSED COMPLETELY. IF THE DOOR HAS BEEN OPENED DURING THE CYCLE, THE STERILIZER WILL DISPLAY: DOOR IS OPEN DURING THE CYCLE! IF SO, PLEASE PRESS START KEY TO CANCEL THE ALARM, THEN CLOSE THE DOOR COMPLETELY AND RESTART.


6.7 Starting a program

Close the door completely, and press START button to start a working cycle, the LCD will display the cycle state. The sterilizer will heat, sterilize and drying the instruments automatically for you. The whole process will take 20-50 minutes. It depends on the object being sterilized, the initial temperature, and the program you selected.

Process Display:

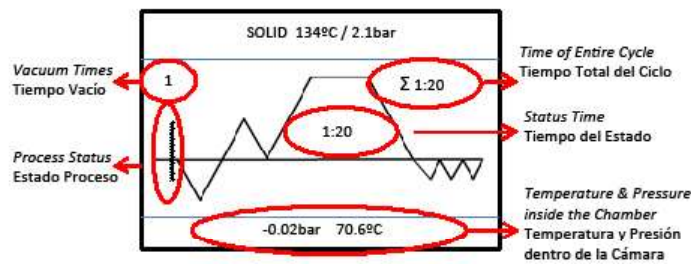Cycle finished:

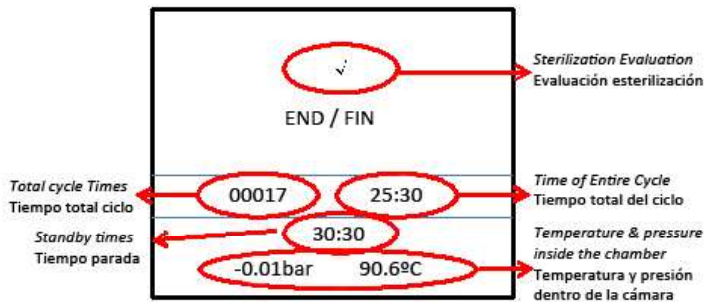CAUTION: DO NOT COVER OR PUT ANYTHING ON THE MACHINE TO KEEP HEAT WELL VENTING.

6.8 End of the sterilization cycle

When working cycle finished, the LCD Screen will give you a sound of alerting, and then you can open the door and take the sterilized instruments out.

WARNING: Do not try to open the door when the pressure value is higher than 0.1 bar.

When the door of sterilizer opened, the program will return to the initial state, heat-preserving and waiting for next sterilizing cycle, before start a new program it will be kept in a heat-preserving condition all the while.

CAUTION: AFTER STERILIZING HAS BEEN FINISHED, PLEASE USE THE HANDLE TO TAKE THE TRAYS OUT FROM STERILIZING CHAMBER. IT WILL BE BETTER TO STORE STERILIZED INSTRUMENT AFTER THEY HAVE BEEN COOLED DOWN TOTALLY.

6.9 Apagado


If you finished the sterilization, please turn off the power switch. The power switch light will be off, and close the door but do not lock it. If you are not going to use it for a long time, please unplug the power cord.

ATTENTION: DURING THE STERILIZATION, WE SUGGEST THAT YOU USE THE INDICATOR TAPE. PUT THEM IN THE CHAMBER IN ORDER TO ENSURE RELIABILITY OF STERILIZATION.

7. Logs Reading software


1. Double click the desktop shortcut (after the program is installed)
2. The Reading software operation screen will appear.

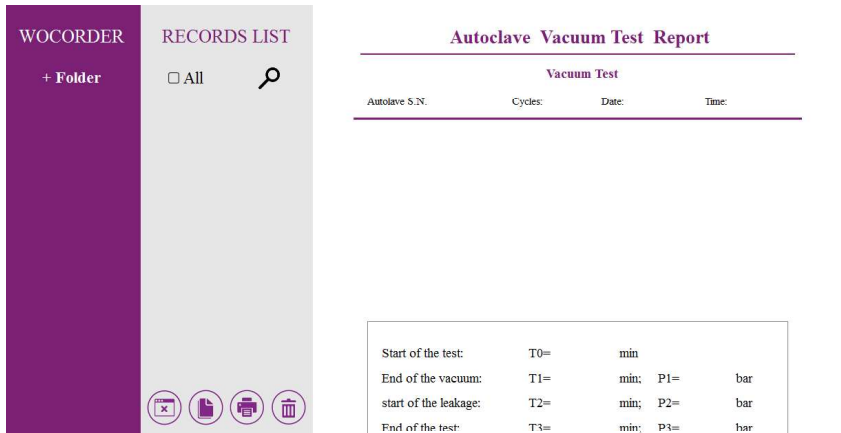3. Click the left mouse button.


4. The following dialog box will appear:

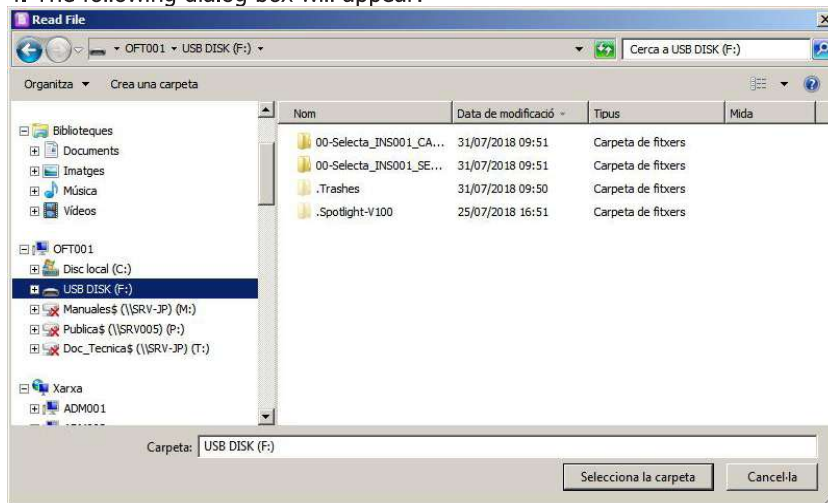
5. Click the **"Choose this folder"** button, then return to the Reading Software main screen. Or click the **"+ Folder"** again, and you can continue choosing another folders.

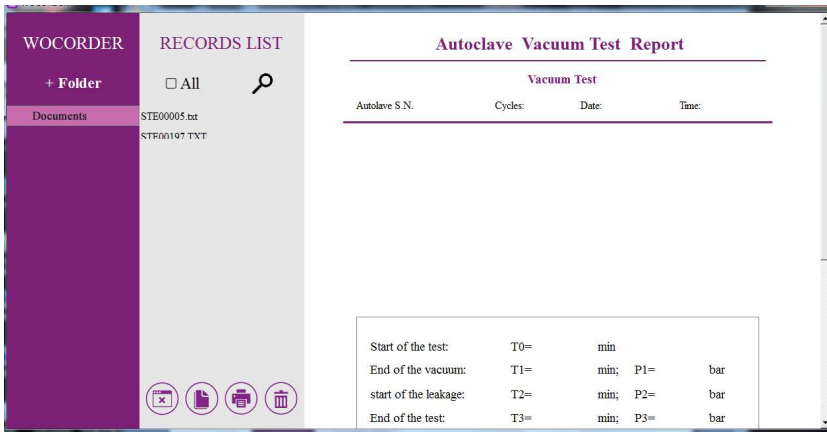
6. Choose one folder in the left column, then it will show you the records in the middle column. If you choose another folder at the same time, they also can be shown in the middle column.

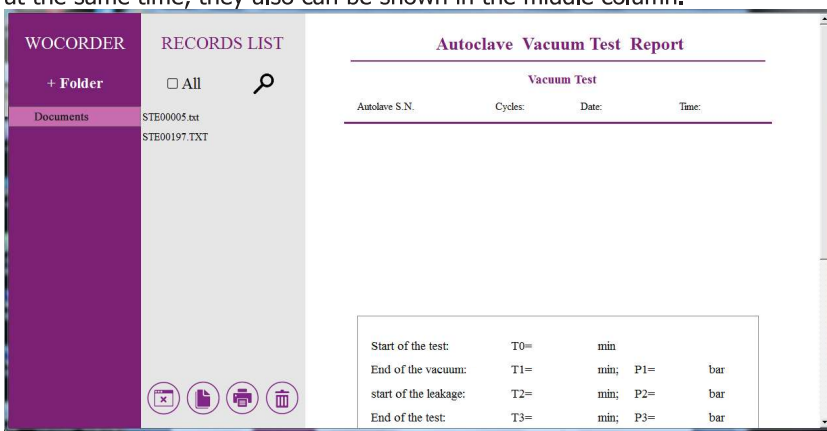
7. When you choose one of these records in the middle column, the log details will be shown in the right screen. If you choose all of the records, it will only show you the first one.

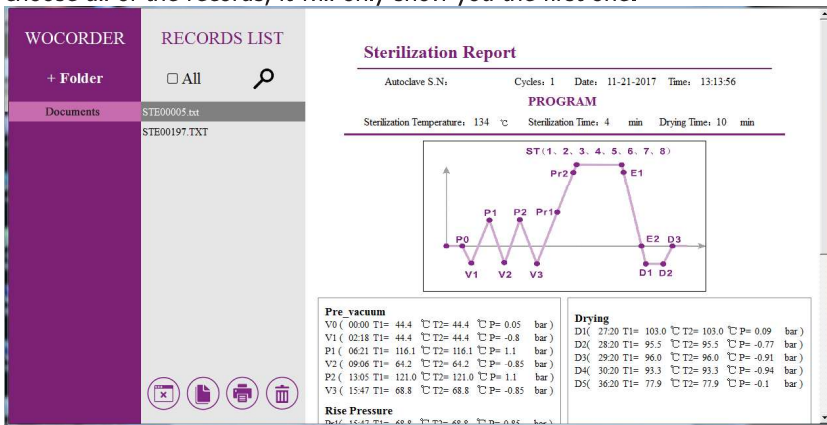
8. Press this button, a dialog box will show up to confirm whether you want to delete all these records in the middle column: choose "yes" and delete all; or choose "no" to return.

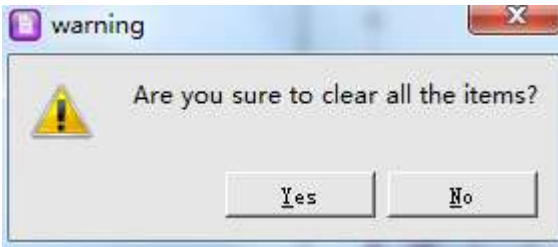
9. Press this button in the middle screen, and copy the chosen records to any document you want.


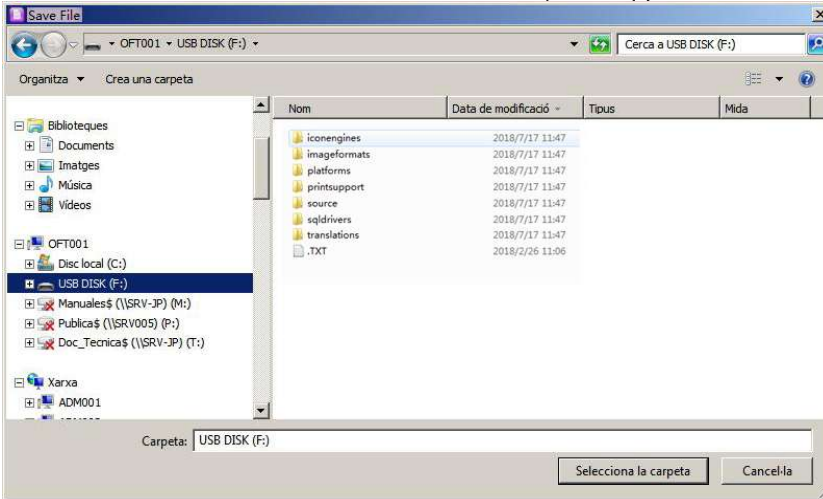
10. Press this button, print the chosen record, a print preview screen will jump out. Press  on the upper right corner, then you can print the record.

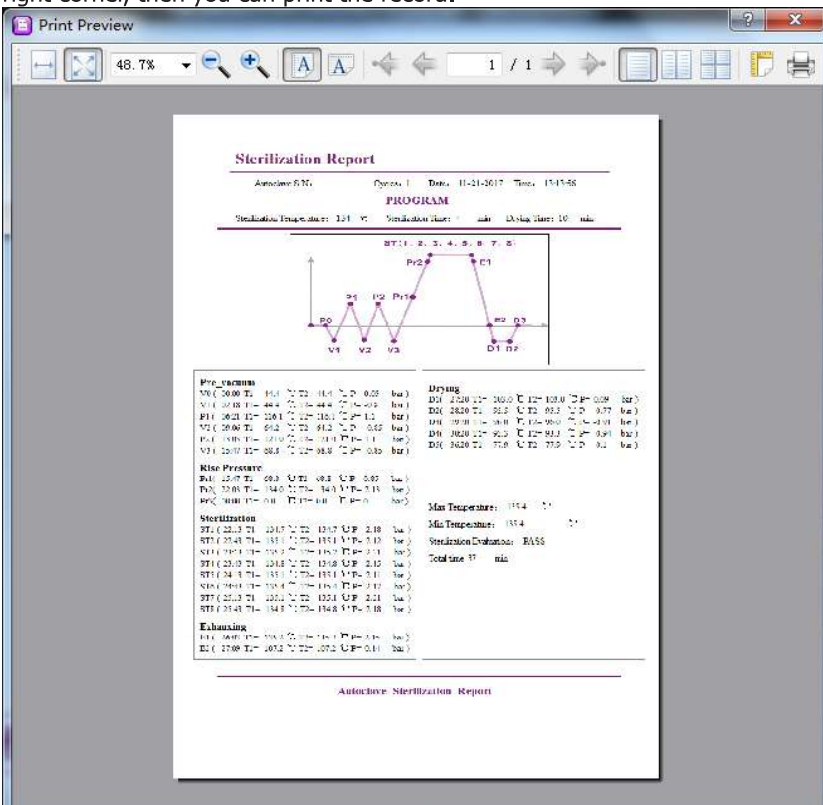11. Press this button, to delete the chosen record.

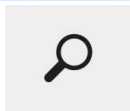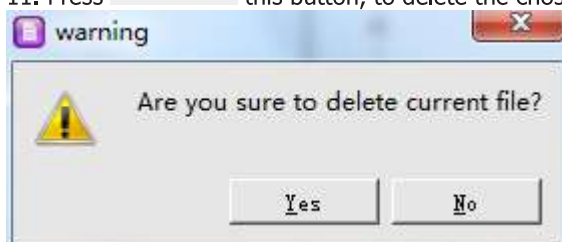12. Press this button, so that users can search the record.

13. The hidden buttons (they appear when the mouse is close) are: the left and right scroll bars in the right screen and the three icons in the upper right corner of the software: minimized, maximized, and closed.

14. The record is shown in the right screen.

8. Maintenance

A general review made by a technician is recommended every two years or every 2500 cycles.

8.1 Maintenance schedule

Daily:

- Clean the door gasket. This action can be performed by the user.
- Clean the sterilization chamber. This action can be performed by the user.

Weekly:

- Clean the chamber, trays and support. This action can be performed by the user.
- Clean the drain filter. This action can be performed by the user.

Monthly:

- Clean the water tank. This action can be performed by the user.

Yearly:

- General revision and maintenance. This action must be carried out by a specialist technician.

When required:

- Replace the door gasket. This action can be performed by the user.

8.2 Daily cleaning

Cleaning the door gasket

The door gasket and the mating surface should be wiped off clean each day with a clean, damp cloth. Do not use abrasive cleaners on the gasket or mating surface. Use warm soapy water for keeping marks of sterilizer persistent, but ensure any soap residues are completely removed by wiping both the gasket and the vessel again with water using a lint free damp cloth.

WARNING

Refer to qualified personnel for servicing: Never use a wire brush, steel wool, abrasive material, or chloride-containing products to clean door and chamber assembly. "Caution hot surface. Avoid contact." ensure that the sterilizer is cooled down fully before cleaning to avoid burns.

Cleaning liquid loads

Cleaning After Liquid Loads Biological media tends to boil at a higher rate than other liquids during venting. This causes media to be spattered inside the chamber. Therefore, the chamber must be cleaned daily when you are sterilizing media. Cleaning as follows:

- Allow unit to be cooled down.
- Wipe out chamber and door with a clean, damp cloth.

8.3 Weekly cleaning

Cleaning Chamber, Trays and Rack

At least once a week, the trays and tray rack should be removed from the sterilizer chamber. The trays, tray rack and chamber should be thoroughly cleaned to remove any deposits from the surfaces. Clean the trays, rack and chamber (especially the bottom of the chamber) with appropriate antibiological cleaners. Wipe all residues from the surfaces with a dampened, lint-free cloth.

WARNING

To prevent from collection of mineral deposits and corrosion of chamber components, use distilled or deionized water only as specified. Clean chamber after each use if sterilizing saline solutions.

Cleaning water draining filter

Water draining filter might have been jammed by some dust because of use for a long-term, so effect of vacuum and drying would be influenced. Some tiny impurity might be deposited on the filter after a long-term use, blocking the filter, so as to influence the effect of the vacuuming and water discharging. The kinds of impurity come from smeary dust on the instruments being sterilized or some calcification in the water. Keep cleaning of the inside chamber in order to make life-time of filter much longer; please take the following advice for consideration:

- Use eligible distilled water;
- The instruments should be cleaned before placing in; it is good to use specified packing for the instruments with oil or other impurity, don't forget to seal up.
- Rotate the water filter which composed by filter net tube(A) and filter holder(B) inside the chamber. Cleaning the part A and B, ensure that there is no any dirty thing on it.(We suggest you to clean by ultrasonic cleaning machine). Then set it back and rotate to the hold bottom of the chamber.

A B

8.4 Limpieza mensual

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There are some impurities and some toxins had been left behind in the reservoir because distilled water stored for a long time, so you need to drain and clean regularly. According to the picture which shows below, loose the screw by screw-driver, and open the cover to clean inside.

ATTENTION

BE SURE THAT USE DISTILLED WATER PROPERLY IN ORDER TO EXTEND THE STERILIZER LIFE-TIME. DO NOT RAVE ABOUT THE STERILIZER WHEN THE TANK HAS BEEN FILLED.

8.5 Other maintenance

Changing Door Gasket

A plain screw driver without sharp head is needed.

Disconnect the sterilizer from the power supply. Ensure that the sterilizer is cool and depressurized.

- 1) Hold verge of the seal by one hand softly, and another hand should be inserted the screwdriver into the gap between gasket and door, take out the seal slowly.
- 2) Once you take out one part of the seal, you can draw out the whole seal slowly. After taking out the seal, please check and clean the groove of gasket, so does the gasket, please replace the gasket if there is some damage.
- 3) Fix the clean gasket in initial door groove. Attention: the gasket should be imbedded in the groove equably. At first, please imbed the 4 spots equably into groove when fix the gasket, and then embed the other parts. After that, press the gasket equably by hand.
- 4) Attention: the inner edge of gasket may be ectropium during embedding it in the door groove, at this time; you'd better to tight it back to the groove by using the screwdriver carefully.

8.6 List of alarm codes

The sterilizer will show an Error information when it has some problem. Here is the Error Alarm List:

CODE No.	Error Information
	Abnormal exit
01	Temperature sensor fault (steam generator)
02	Temperature sensor fault (out chamber)
03	Temperature sensor fault (inner chamber)
04	Sterilization fail
05	Pressure not exhausted
06	Door open in the cycle
07	Operation over time
08	Over pressure
09	Over temperature
10	Pressure and temperature do not match
11	Vacuum process fail
98	Power fail

You have to press Start/Stop button for more than 3 seconds to exit the error window. If the error occurs again, contact you service.

8.7 Items to inspect during annual maintenance

1	Exterior	The exterior of sterilizer should be tidy and mustn't have the disfigurements, such as deflection, hollowness, collision, nick, sharp edge, and so on.
2	Cover plate	The cover plate should be assured to disassemble easily in order to repair the equipment.
3	Digit and letter	The digit and letter in screen should be legible.
4	Electroplate components	Revise status & operation
5	Printer components	Revise status, operation and printing quality
6	Door safe lock	On the normal condition, if the sterilizer door hasn't been locked tightly, the program can not start.
7	Chamber Pressure	The sterilizer should ensure that the door can't be opened when chamber is pressurized
8	Safety Valve	The sterilizer must install a safety valve, safety valve opening pressure 2.5 barg, and automatically open when reaching the set value.
9	Sterilizing Program	Sterilizer should have the pre-established program about 121°C and 135°C, dressing and instruments.
10	Controlling system	The control system in sterilizer should limit the steam which in the chamber be controlled at the highest average temperature in $\pm 3^\circ\text{C}$ of pre-establish station. And ensure the temperature value accords with the pressure controlling value.
11	Timing control	Able to timing control of the sterilization and drying, and the error should not be greater than 10% of the preset value.
12	Button and Switch	Buttons and switches should be flexible and reliable on the sterilizer.
13	Indicator and display	The indicators and displays of sterilizer should show the states of every sterilizing procedure exactly. Under the normal situation, sterilizer should indicate: a) Chamber temperature b) Chamber pressure c) Sterilizer working state d) Water level state e) State of door

14	Quantum of leakage	In the condition of the vacuum, the sterilizer shouldn't leak 0.0013Mpa within ten min.
15	Leakage forbidden	The sterilizer can't leak under the work pressure
16	Protective earthing impedance	The impedance between protective earthing point of the power input faucet and protective earthing can be touched all metallic parts, doesn't over than 0.1Ω.
17	Empty load	For all loads except hollow load A, the presence of saturated steam in the usable space and the load is deemed to have been achieved when, throughout the holding time, all temperatures measured in the usable space and the load: (Attention: the theory of steam temperature is accounted by measuring pressure, which can be considered the test temperature.)
		are not lower than the sterilization temperature;
		are not more than 4 K above the sterilization temperature;
		do not differ from each other by more than 2 K.
		The usable place temperature during the no-load can not over than the scope of highest temperature.
18	Hollow load	For hollow load A and B, in order to confirm the presence or absence of saturated steam, discriminate whether the indication system change in accordance with the system manufacturer predetermined color.
19	Dryness, solid and wrapped load	For wrapped load, any remaining moisture should not lead to wet packages and shall not result in detrimental effects on the sterilizer load. Any remaining water droplets on the inner side of pouch should evaporate within 5 min.
		For solid load the moisture content should not exceed 0,2 %

9. Transportation and storage

9.1 Preparation before Transportation and Storage

Shut off the power switch, unplug the cord, and make the sterilizer been cooled down completely.

9.2 Draining

Drain water from reservoir and the condensate collector completely: insert the joint end of the attached tube to drain connection. (The spout on the left is the water spout used for the "used-water out", the one on the right is used for the "clean-water" drain spout).

Used-Water out Clean-Water out

9.3 Conditions For Transportation And Storage

Temperature: -5°C ~ +55°C

Relative Humidity: ≤85%

Atmospheric pressure: 500HPa~1060HPa

10. Warranty

This product is guaranteed for one year. The warranty does not cover damages caused by improper use or by causes beyond J.P. SELECTA, s.a.u.

Any manipulation of the equipment by personnel not authorized by J.P. SELECTA, s.a.u. automatically voids warranty benefits.