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Ejemplo:

```
R,Biachi,, 0, 43, 1, TP, 136.00, g/L, 0794
R,Rossi,,1, 40, 2, TP, 147.00, g/L, 0765
R,Tony,,0,6,3,TP,35.00,g/L,069d
TRANSFER FINISH
```

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1 Packing list and accessories


1.1 Paking list

1. Analyzer
2. IEC Power cord
3. Anti-dust cover.
4. 2 units of 58 mm wide thermal paper roll (Ø 50mm)
5. Spare lamp and lamp holder
6. Set of spares for tubing
7. Serial 9 pin communication wire.
8. 5x20mm fuse rated T3.15 250V
9. Serial mouse.
10. Instruction manual.

1.2 Accessories

58 mm wide thermal paper roll.

4401156

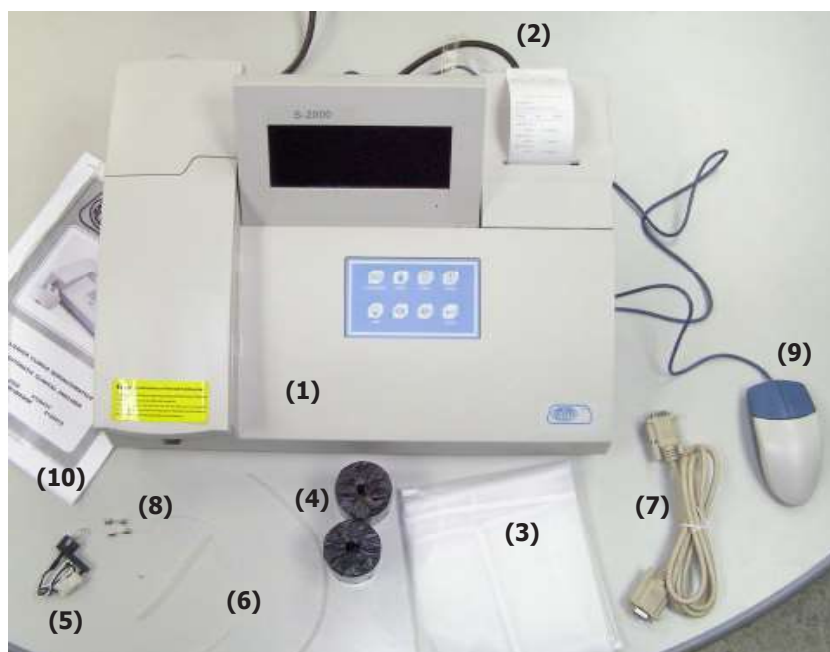2 Analyzer overview

The S-2000 is a photometric equipment designed, mainly, for clinical analysis uses.

The control software (Windows CE based) brings an easy and intuitive usage.

In order to allow a fast start-up, 47 most popular tests are preprogrammed on the unit. The user can add up to 130 tests. Test modes are:

- End point
- Bichromatic
- Lineal regression
- Multi-standard calibration
- Two-point
- Kinetics

The S-2000 needs a mouse (included) or a keyboard (not included) to be operated.

"Quick keys" allow a fast access to most common functions.

3 Safety cautions and potential hazards

3.1 Symbols used on this manual

Important information

- Information or tips to obtain good results or for an optimal operation of the equipment.
- Information or tips to extend the life of the equipment or to avoid the degradation of some of its components.

Warning

Warning of a dangerous operation. Follow indicated instruction.

Hazard of electrical shock

Electric hazard for entering in areas marked with this signal, or for following the manual instruction marked by this icon.

Follow the instructions for the described operation.

Hazard of burns to be in contact with hot surfaces

Temperature in this area can exceed 60°C. Use heat protected gloves for the described action.

Follow the instructions for the described operation.

Hazard of contact with corrosive substances

Hazard of contact with sodium hydroxide and strong acid, substances that attack skin, clothes and other materials easily.

Follow the instructions for the described operation.

Biohazard

Possible contact with infectious or biological dangerous material.

3.2 General hazards

Using S-2000 needs some cautions:

• Operator Qualification

Only medical qualified personnel must use and take conclusions from the results given by the S-2000

• Service Technician Qualification

Contact with your dealer to get qualified service for this unit.

• Electrical Hazard

- Do not use the unit without the external covers.
- Unplug from mains to run any maintenance or service operation.
- Mains voltage on the main electronic boards.

• Lamp area

Lamp becomes extremely hot during operation. Never touch the lamp when it is on!. Avoid to look directly into the light. It could cause eyes damage.

• Chemicals and reagent

Working with S-2000 needs chemical manipulation. Get Chemicals manufacturer information for safety operation. Avoid direct body-contact with reagents or chemicals.

Biohazardous Materials

Treat all **samples** and **waste** as potentially biohazardous and infectious.

Avoid direct contact with samples, waste solution and/or solid waste.

Both should be handled as potentially biohazardous.

Use gloves, glasses lab coat,.. as a general rule to operate S-2000 with real samples.

If any sample is spit on the instrument, use qualified procedures for decontamination and disinfection.

Avoid direct contact with waste solution and/or solid waste.

Both should be handled as potentially biohazardous.

Biohazardous parts

Avoid direct contact with the sipper tubing and all parts of the sample flow path. Treat these areas as potentially biohazardous and /or infectious

3.3 Additional Precautions

• Flammables

Avoid using dangerous flammable material around the instrument.

• Accuracy/Precision of the Measured Results

For proper use of the instrument, measure control samples and monitor the instrument during the operation.

An incorrectly measured result may lead to an error in diagnosis, thereby posing a danger to the patient.

Treat all reagents according to the manufacturer's recommendations. Refer to the reagent kit box and package inserts, or product information sheets for specific instructions.

Make sure that the sample/reagent mixture does not contain any blood clots, dust or other insoluble contaminants. If insoluble contaminants are contained in the sample, correct measuring values may not be obtained.

• Application

The instrument is designed for clinical chemistry test analysis using water-soluble samples and reagents.

Please note that other types of analysis may not be applicable to the instrument.

• Operation and Maintenance

During operation and maintenance of the instrument, proceed according to the instructors and do not touch any parts of the instrument other than those specified.

Never leave a Reagents/sample mixture in the flow cell for longer than necessary. Always clean the flow cell after a batch of measurement and keep the flow cell filled with distilled water when not in use.

Verify the front covers closed while the instrument in operation.

Avoid touching the mechanism, such as the sipper mechanism inside the instrument, while the instrument is operating. This may cause operation stop or damage the instrument.

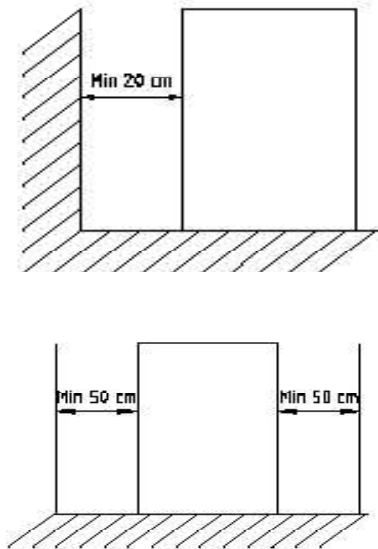
Fig 4.1 and 4.2 let free spaces around the S-2000.

4 Installation

Unpack the instrument, removing it from its plastic bag. Keep the boxes and packaging material until being assured that nothing is missing, that the unit has no (internal) transport damage and are working properly.

4.1 Instruments location

Choice a properly location for the S-2000:

- Flat surface, levelled, stable and suitable for its weight.
- Avoid direct sun light and dusted ambient.
- Far from equipment that can produce vibration. Eg: centrifuges, etc...
- Not far (1.5m) from a power socket.
- It is suggested to have a waste near the S-2000.
- Keep free 20cm on the rear part. See figures 4.1 and 4.2
- To operate comfortable, let free 50cm on both sides. See figure 4.2
- Working ambient: Temperature: 15-32 °C, humidity: 20% - 85%

4.2 Connect the instrument to AC power

Plug the power cord into the instrument. Plug the power cord into an AC electrical wall outlet.

AC110V ~AC250V 50~60 Hz 150W max

ATTENTION! Cautions

AC power outlet must be connected to ground.
Shut (Power) off the power if smog, special smell or special sounds from the inside of the instrument, and contact your supplier.

Hold the connector when you pull off the power cord. Do not pull the cord directly.

Figure 4.3 Adjust display position.

Figure 4.4 Removing protective covers

4.3 Paper Installation

- Remove the paper cover by pulling up at the top end of the cover.
- Unroll about 10 inches of paper and place the roll on the table behind the instrument.
- Feed a cleanly cut edge of the paper from the back into the printer.
- A ragged edge or wrinkled paper will be difficult to load and could cause a paper jam.
- Feed just over 1 inch of paper in and then press the feed paper key several times to feed over 3 inches of paper. Drop the roll of paper into the printer paper well in back.
- Feed the paper through the paper slot of the cover, then close the cover.

4.4 Flow cell and cuvette system

- Remove flow cell cover
- Insert the flow cell or cuvette into the flow cell holder
- Make sure the surface with characters of the flow cell is faced to user.

Note!


Keeping the surface of the flow cell or cuvette clean. Don't touch it by fingers.

5 Unit description


5.1 Front View

- | | | |
|-----------------------|---------------------|----------------|
| 1. Pump cover:* | 2. Flow cell cover* | 3. Sample tube |
| 4. Sample button | 5. LCD Screen | 6. Hot-key |
| 7. Power-on indicator | 8. Printer cover | |

* Open the cover, just to change the lamp or maintain the pump.

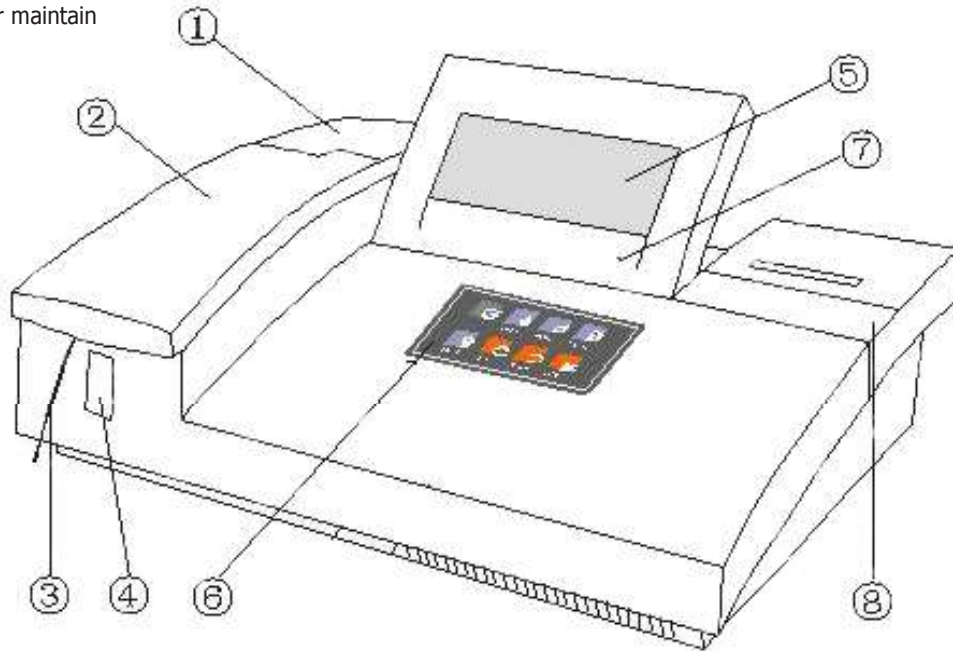5.2 Quick keyboard

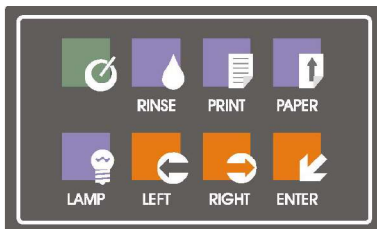
5.1 Quick keyboard

- Soft keyboard:** Open or close soft keyboard.
- Rinse On- Off:** Pressing this control will switch on the sipper pump in a continuous mode. Pressing the key again will stop the pumping action at the end of a complete cycle.
- Feed Paper:** Control the printer paper for the embedded printer.
- Print:** Set the printer enable or disable.
- Lamp On-Off:** Turn off or turn on the lamp.
- Left:** This key is used to go backwards in a procedure.
- Right:** This key will advance the cursor to the next position.
- Enter:** This key is used to confirm an entry via the keyboard.

5.3 Rear view

- | | | |
|---------------------|----------------------|-------------------|
| 1. Power switch | 2. Fuse | 3. AC Mains port. |
| 4. Fan | 5. Printer interface | 6. Mouse |
| 7. RS-232 interface | 9. ext. keyboard | 10. Not Used |
| 11. Drain Port | 12. Pump | |

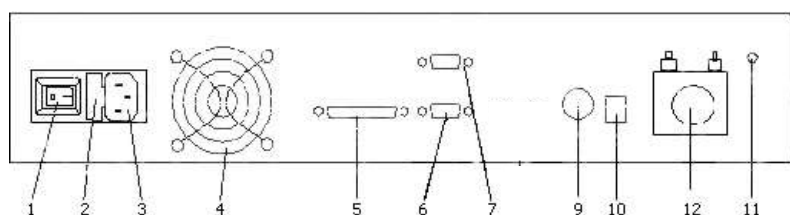5.4 Flow cell and cuvette

The following sketches will help you locate and identify the flow cell parts.

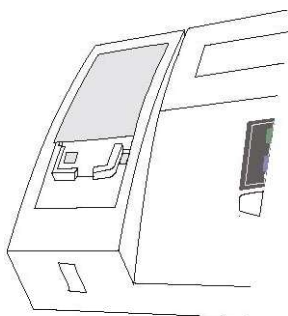
Figure 5.2 Accessing to lamp, flow cell and sipper pump.

Sample. Flow cell. Sipper pump.

5.5 Lamp and sipper pump

1. Remove the lamp cover and flow cell cover.

Figure 5.3 Sipper pump

CAUTION

The source lamp becomes extremely hot during operation;
never touch the lamp when it is on!

Look directly into the light path of the lamp may cause eyes damage.

If the lamp needs to be changed, always switch off the lamp by switching off the instrument and then wait until the lamp has cooled down.

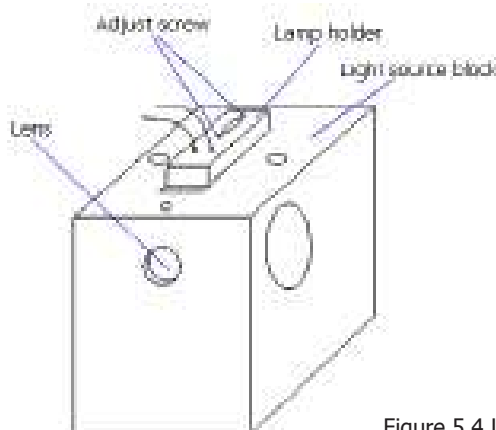
Figure 5.4 Lamp holder

6 Operation

6.1 Power on

Turn the instrument on using the power switch on the rear panel. The display takes 10s to appears.

Fig 6.1

The instrument will be waiting the lamp stable in 2 minutes. Then the Fig 6.1 will display figure:

Raise a test tube of distilled water up to the aspiration tube.

Press "Rinse" key to rinse the flow cell for about 20 seconds, then press the "Sample Button" to aspirate the distilled water.

The S-2000 will now read and save the absorbance reading for water for all filters in the machine.

If the absorbance reading is too high or too low, the Fig 6.2 will be displayed.

Screen 6.2

Select "No", to ignore the deviation above and move to test menu.

Select "Yes", to repeat the measure.

6.2 Basic operating

Mouse

The instrument uses normal mouse for PC (RS-232 interface)

Keyboard (option)

Connect the ext. Keyboard to the instrument when the analyser powered off. You can use the keyboard just like to operating a PC keyboard.

Soft-keyboard

Use the hot-key or click the icon (soft-key) to pop a soft-keyboard, (screen 2.5 Soft keyboard) and click the icon again, the soft-keyboard will be closed.

Soft-key

Icon (soft-key)

Click here when dragging

Screen 6.3 Soft keyboard

The soft keyboard is used for entering data which needs characters or numbers.

If you want to move the keyboard, please click its upper blank area and drag it away.


6.3 Main menu and system parameters setting

After the instrument initialization, the screen 6.4 is displayed.


Select "Main" to set the main parameters for the unit.

Screen 6.4 Fixed test menu

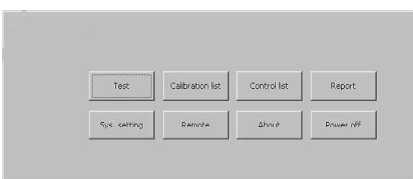
Screen 6.5 Main menu

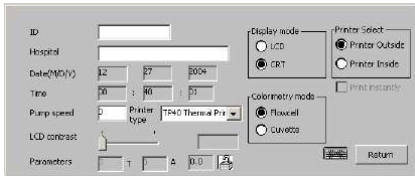Screen 6.6 System settings menu

System settings menu:

1. Setting time and date
2. Select pump speed (Factory set at 6)
3. Adjusting contrast of LCD
4. Printer Select
Printer Inside: Thermal printer is selected.
Printer Outside: External printer.
5. Enter the Hospital or laboratory name

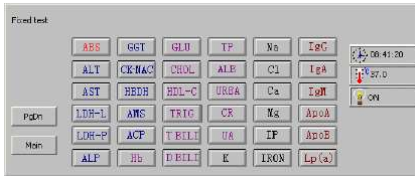Screen 6.7 Fixed Test menu

6.4 Programming test parameters

S-2000 Analyser has a built-in menu of 130 assays. 47 assays of most commonly applied data reduction formulae are preprogrammed.


6.4.1 Modifying stored programs

Click "Test" button in main menu. The screen 6.7 will be displayed.

1. Click the button which will be modified, such as "UA", The screen 6.8 will appear.

Select "Edit" to modify the parameter of the test. Screen 6.9.


2. Modifying parameters according to the reagent instruction.
3. Click "Save" button to save the modification, or click "Cancel" to ignore.

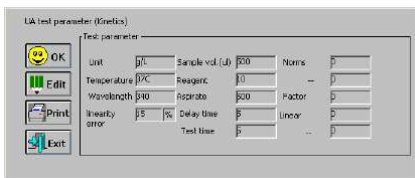Screen 6.8 Test Enter menu

6.4.2 Programming a new assays

1. Select a free button of assays by clicking "Pgdn" button, then click the button.
2. Input parameters
3. Save the input

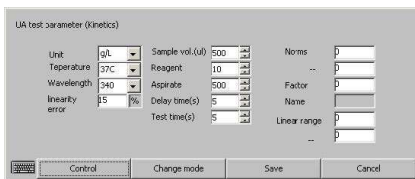Screen 6.9 Programming menu

6.4.3 Calibration parameter

Select "Calibration", "Calibration setting" window will be open.


Screen 2.10 Calibration setting menu

1) Calibration for linear tests


- Standard number means that how many standards will be measured.
- The measured results will be averaged. Input the number by clicking the window.
- Repeat number: This parameter determines how many standards (with the same concentration) will be used for calibration. The results of multiple measurements will be averaged by the system.
- Select "OK" to confirm the selection.
- Select "Exit" to go back to above menu with modify

2) Calibration for non linear tests

- Standard number means that how many standards will be measured.
- The standards are different concentrations. The measured results will be a calibration curve. Input the number by clicking the window.
- Repeat number: This parameter determines how many standards (with the same concentration) will be used for calibration. Possible inputs are 1, 2 and 3. The results of multiple measurements will be averaged by the system.
- Select "OK" to confirm the selecting.
- Select "Exit" to go back to above menu.

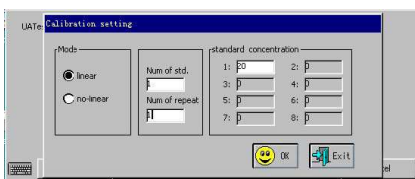Screen 6.10 Calibration menu.

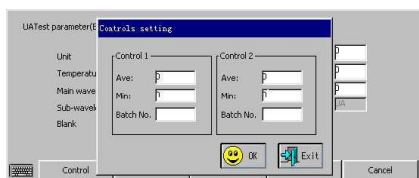
Screen 6.11 Controls setting menu

6.4.4 Control parameter

Select "Control" for programming control. (Screen 6.11)

- | | |
|--|--------------------------------------|
| a) Control 1 and control 2: | Two controls could be defined. |
| b) Mean value: | Input the mean value of the control. |
| c) Minimum: | Input the minimum of the control. |
| d) Batch number: | Enter the control batch number. |
| e) Select "OK" to confirm the programming or Select "Exit" to go back the above menu without modifies. | |

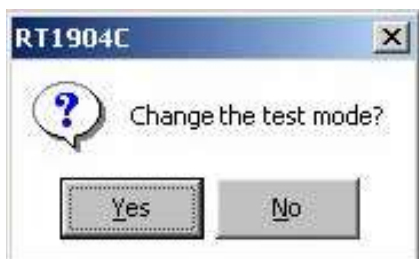6.4.5 Selecting a Mode

Into the menu of the "Programming" (screen 6.9), click "Change mode". And then select "Yes", the screen 6.12 appears:

a) Endpoint

After a required incubation time, the reaction reaches its maximum color development, or its Endpoint. The color of the reaction will remain stable and constant for certain period of time. During this time a measurement in Absorbance is taken.

b) Endpoint Bichromatic

Bichromatic is actually endpoint determination measured at two different wavelengths. The advantage is that some unwanted interference due to the color of the serum sample can be eliminated by subtracting the measurement done with the second wavelength filters. This can be done when you select a sub-wavelength in "Parameter change" menu. It was described in 6.4.1.

c) Endpoint determination against a standard Curve or against a Multi-standard Curve


This is a normal endpoint determination, in which linearity may not be perfect throughout the whole range that is interest. As a result it is calibration curve (a number of standards) instead of a single standard.

d) Kinetic

Kinetic determination are tests whose enzyme activity is measured continuously over a certain period of time and finally referred to 1 minute.

e) Two point

After a pre-defined incubation time, two absorbance readings are taken at different and programmed time interval. From the delta absorbance between the two measurements the final result is calculated. Each of the measurements taken is a mean of a large number of readings taken in quick succession.

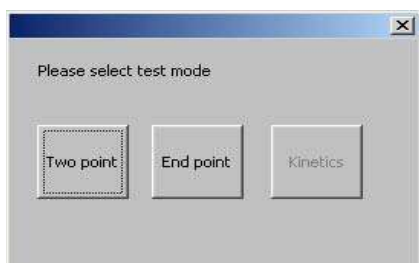Screen 6.12. Change mode menu


6.4.6 Save the Change

After select the measure mode, Select "Save change" to save the selection. The menu will change to the "test enter" menu. Select "Enter" to perform the test procedure.

7 Performing tests

Select the test to be performed in the menu of "test overview", and select "OK" to go to the menu of "test procedure" in the menu of "Test enter".

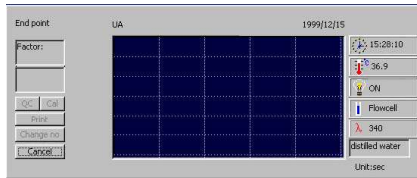
Screen 7.1 Test procedure.

7.1 Test procedure menu

The "test procedure" menu is displayed as follows:

- | | |
|-------------------------|---|
| 1. Factor window: | Factor |
| 2. Result window: | Test result |
| 3. Control: | Perform the control test. |
| 4. Calibration: | Perform the calibration test. |
| 5. Print: | Print the test result. |
| 6. Change No.: | Change the sequence number tested. |
| 7. Time: | Display the time. |
| 8. Temperature: | Display the temperature of the flow cell. |
| 9. Lamp: | Show the status of the lamp: On/Off. |
| 10. Wavelength: | Show the wavelength selected. |
| 11. Instruction window: | Display the instruction information. |
| 12. Return: | Return back to above menu. |
| 13. Reaction window: | A curve will be displayed to represent the reaction |

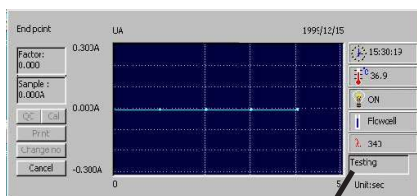Screen 7.2

Instruction window

7.2 Running a test

To do a test, follow the instructions shown on the bottom right corner "instruction window".

1. Instruction window: Aspirating distilled water

Put a tube with distilled water under the sample tube and press sample key (figure 7.3). The measure result will be displayed on the "result window" and on the "reaction window".

2. Instruction window: Aspirating reagent blank or aspirating sample blank


Hold the reagent blank or sample blank under the sample tube and press the Enter key or the sample button. The programmed solution will be aspirated and measured.

3. Instruction window: Aspirating Sample

Hold the sample under the sample tube and press the sample button. The programmed solution will be aspirated and measured.

4. Calibration: It will be described in Chapter 8

5. Control: It will be described in Chapter 9

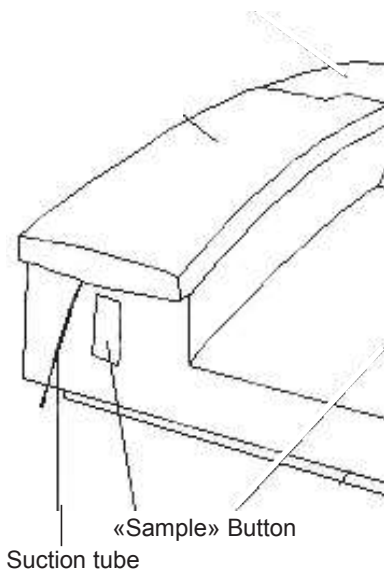Suction tube

On pressing «Sample» button, the sample is suctioned.

7.3 Changing test

When changing over from one test to another, washing the flow cell is required to prevent reagent carry over, causing wrong calibration and measurement result.

To wash the flow cell:

- Place a container of distilled water under the sample tube. Press the **Rinse** key. The aspiration/flow cell system will be washed continuously now, until the **Rinse** key is pressed again. Let the flushing continue for about 30 seconds.

Fig. 7.3. Sample entry for its analysis.

8 Calibration by standard

8.1 General

Certain tests require calibration by 1 or more standards. Calibration results are stored in the flash memory of the analyser and can be used again for future measurement. Calibration can be programmed and subsequently done for different measurement mode: Endpoint or Twopoint. For Endpoint, the calibration could be with/without Reagent Blank and/or Sample Blank. For Twopoint, the calibration could only be with or without Reagent Blank.

The calibration by standard should be performed depending on the stability of the reagents, batch of the reagent, and type of test. The new ones will overwrite the previous calibration data in the flash memory.

8.2 Program calibration

The calibration parameters could be programmed in the menu "Parameter modify", see the screen 8.1. The setting was described in paragraph 6.10

8.3 Perform linear tests

For linear tests, the relation between Absorbance and Concentration is a (linear) straight line. This test usually is calibrated with one or more standard(s) of 1 concentration and with or without Reagent Blank.

Select "Calibration" in the "test procedure" menu:

1. Measure the standard

Instruction window will appear: Standard "#"


Hold the standard sample "#" under the sample tube and press Sample button. The programmed volume will be aspirated and measured. The process and result of the measurement will be displayed.

2. Repeat the measure

Instruction window: Standard "#", Repeat #

3. Abs / Concentration

The result will be displayed on the result window, and the standard curve will be displayed in same time. User can save and print the result.

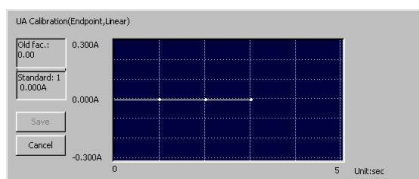
Screen 8.1

8.4 Perform the calibration of non linear tests

Non linear tests are those tests where the relation between Absorbance and Concentration is not linear. These tests are usually calibrated with more than 1 standard of different concentrations. The analyser has a special arithmetic to link the calibration points to a calibration curve. The measurements for patient or controls are then interpolated automatically depending on the curve.

The calibration can be run similar as 8.3

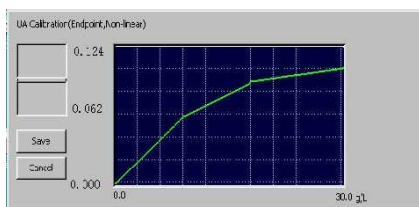Screen 8.2

Screen 8.5 "List screen"

8.5 List calibration

Select "Calibration" in the "Main menu", to list the calibrations. (Screen 8.5)


- 1) Select the test which one will be checked. Click "OK".
- 2) The result of the test will be displayed on the screen.
- 3) Print the result by selecting the "Print".
- 4) Select "Exit" to go back to the above menu.

9 Control

9.1 General

The control program gives the possibility of storage and statistical evaluation of controls measured on the instrument.

For all 130 tests programmed in the unit, 2 controls can be defined. The result of controls for 130 tests of 365 days could be stored automatically.

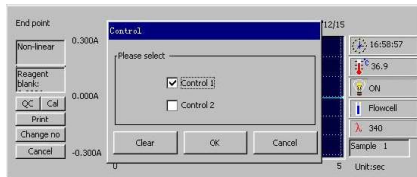Screen 9.1

9.2 Program controls

The control parameters could be programmed in the menu "parameter modify", see the screen 2.15. The setting was described in paragraph 6.4.4

9.3 Measurement of controls

A control can be run at any time during a test run, after the test has been calibrated.


Select "Controls" in the "test procedure" menu, a window will be open as following: (Screen 9.1)

1. Choice "1" if only measuring 1 control, and choice "2" if 2 controls will be measured. Select "Enter" to confirm the choice.

2. Instruction window will display: Control "#"

Hold the control requested by analyser under the sample tube and press Sample button. The programmed volume will be aspirated and measured. The process and result of the measurement will be displayed.

3. The result will be added to the statistical memories automatically. The evaluation of control results can be accessed through the main menu and described in 4.4.

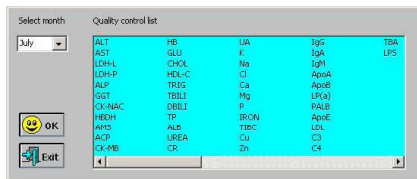
Screen 9.2

9.4 List controls

Select "Control" in main menu, to list the statistics of controls. The list screen (Fig 9.2) will be displayed:


1. Select the test, and the month will be checked. Click "Enter".
2. The statistic of the controls will be displayed on the screen 9.3
3. Print the result by selecting the "Print".
4. Select "Exit" to go back to the above menu.

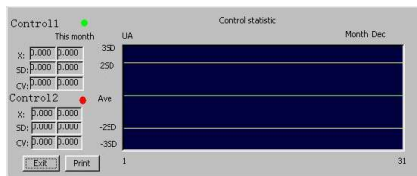Screen 9.3

10 General Specifications

Weight:	10 kg
Dimension:	450mm(L)x330mm(W)x140mm(H)
Power supply:	90-250VAC, 50/60Hz
Power consumption:	150 W max
Fuses:	3.15A, slow-blow type 5x20mm
Operating condition:	15°C to 32°C; RH 85% max
Storage condition:	-20°C to 50°C; RH 85% max
Display:	7" color LCD (640x240, 256 colors)
Measurement range:	0.000-2.500 abs
Resolving power:	0.001 abs (display) 0.0001 abs (calculation)
Interference filters: (nm)	340, 405, 500, 546, 578, 620, 670
Half band width:	<10nm
Lamp:	Tungsten halogen lamp 6V/10W
Temperature control: (Peltier)	25, 30, 37 ±0.1°C and ambient temperature
Flow cell:	30ul quartz flow cell.

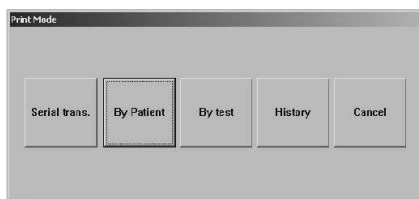Screen 11.1

11 Report

11.1 Summary

Four reports styles are provided:

- Serial trans.
- By patient.
- By test.
- History.

11.2 Serial transfer

Test result can be transferred to PC Computer through standard RS232 using the provided wire:

The serial port of PC should be set as following:
19200 baud, Even Parity, 8 bit of data and 1 bit of stop.

See details on Appendix II:

11.3 Print by patient

Select "By patient", a new window will be open.

1. Update patient information

Select one patient, Line is highlighted. Input name, department, age and sex, click "Update" to validate patient information.

2. Add a patient

New patients can be added, just click "Add" after the patient information input, and a new patient number is attached to it. If you want to test the patient's sample, just change the present sample No. into the patient No.

3. Select print mode

In this part, seven print modes are provided, "General", "Liver function", "Kidney function", "CK enzymes", "Lipids", "Iron" and "Electrolytes".

General: Print all items tested in the test order.

Liver function: The following 18 items are included.

- | | | | | |
|-----------|-----------|---------|-------------|-------------|
| 1. ALB | 2. TP | 3. GLB | 4. ALB/GLB | 5. T BILI |
| 6. D BILI | 7. I BILI | 8. AST | 9. ALT | 10. AST/ALT |
| 11. GGT | 12. ALP | 13. LDH | 14. LDH/AST | 15. GGT/AST |
| 16. TBA | 17. PALB | 18. GLU | | |

Kidney function: If this mode is selected, 6 items are included.

- | | | | | |
|---------|-------|-------|-------|--------|
| 1. UREA | 2. CR | 3. UA | 4. TP | 5. ALB |
| 6. GLU | | | | |

CK enzymes: If this mode is selected, 4 items are included.

- | | | | |
|-------|----------|--------|--------|
| 1. CK | 2. CK-MB | 3. AST | 4. LDH |
|-------|----------|--------|--------|

Lipids: If this mode is selected, 8 items are included.

- | | | | | |
|---------|--------------|----------|----------|---------|
| 1. CHOL | 2. TRIG | 3. HDL-C | 4. HDL-C | 5. ApoA |
| 6. ApoB | 7. ApoA/ApoB | 8. ApoE | | |

Iron: If this mode is selected, 2 items are included.

- | | |
|---------|---------|
| 1. IRON | 2. TIBC |
|---------|---------|

Electrolytes: If this mode is selected, 11 items are included.

- | | | | | |
|---------|-------|-------|--------|---------|
| 1. K | 2. Na | 3. Cl | 4. Ca | 5. Mg |
| 6. P | 7. Cu | 8. Zn | 9. GLU | 10. AMS |
| 11. LPS | | | | |


Input patient NO which need to be printed, patient NO can be input as following:

- Defined numbers. Eg: 3, 6
- Range of numbers. Eg: 1-5
- Various ranges Eg: 1-4, 5-6
- Range and pages Eg: 3, 4, 5-7
- All patients Select the check box "Select all"

Click "Print", to print the report. If some item in the mode hasn't been tested, there will have no correspondent data in the report.

4. Input no-chemistry item

Setting: Directly inputting a no-chemistry item and result is also permitted, this item will be printed in the "General" mode report.

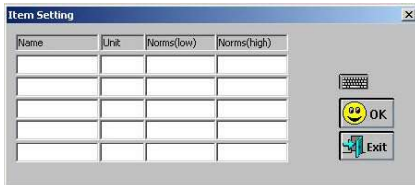Screen 11.2

In the window of "By patient", click "Setting", a window will be open as Screen 11.2. "Test name", "Unit", "Normal (low)" and "Normal (high)" are all need to be input, and each can't be longer than eight characters. All the no-chemistry items will be saved if the Chemistry Analyser is closed correctly.

11.4 Print by test

In the state of Screen 11.1, click "By test", a new window will be open. Select the Item to be printed, then click "Print", all test result of selected item will be printed. If you want print all items, just select the "Select all" check box.

11.5 Print history records

In the state of Screen 5.1, click "History", a new window will be open. S-2000 can store 500 pieces of records at most. In the new window, drag the scroll bar, you can see about all test results of history.

Input the order No. of patient, then click "print", all selected records will be printed out in a history records list. Patient No. can input as following:

- Defined numbers. Eg: 3, 6
- Range of numbers. Eg: 1-5
- Various ranges Eg: 1-4, 5-6
- Range and pages Eg: 3, 4, 5-7

12 Close system

At the end of the working, you must perform the Closing procedure by selecting "Power off" in the "Main menu".

1. The system will store the parameters and the results of tests.
2. The screen will display "aspirating distilled water".
3. Place a container of distilled water under the tube.
Press **Sample button**, Rinse will take few minutes
4. The screen will display "Turn off the System"
5. Turn off the power.

ATTENTION! Cautions

Never power off the system without performing the closing procedure. It will cause that lost of the data and influence the measurement.

13 Maintenance

13.1 General

The S-2000 is a clinical instrument, which requires only a minimum amount of maintenance.

Some general maintenance is necessary to keep the instrument in an optimal condition.

13.2 Cleaning the outside

Always keep the outside of the instrument clean and free dust. You can wipe the outside with a damp (not wet) cloth and a mild detergent.

Wipe the display screen with a soft, non abrasive cloth.

ATTENTION! Cautions

Do not use any type solvent, oil, grease, silicone spray, or lubrication on the instrument.

Place under the «suction tube» a water or cleaning solution to clean the flow cell.

13.3 Clean procedure for the flow cell

The inside of the flow cell should be kept as clean as possible, to assure good and reliable measurements.

To keep the inside clean, follow the following procedures:

13.3.1 Washing when Changing tests

When changing over one test to another, some washing of the flow cell is required. To wash the flow cell:

Place a container of distilled water under the tube. Press Rinse key. The flow cell system will be washed continuously now, until the Rinse key is pressed again. Let the flushing continue for about 2 minutes.

13.3.2 Washing twice a day

The following procedure is needed twice a day, to wash out those residues from the wall of the flow cell, which cause bad filling of the cell:

With the solution of neutral detergent. Let the neutral detergent stay in the flow cell for 2-3 minutes. The distilled water added some detergent is better. The ratio of the detergent is depending on the detergent used. Normally a 1-5% solution is adequate.

Rinse with distilled water without detergent for about 1 minute.

13.3.3 Closing the system by the closing procedure

At the end of the working, you must perform the Closing procedure by selecting "Closing system".

- 1) The system will store the parameters and the results of tests.
- 2) The screen will display " distilled water".
- 3) Place a container of distilled water under the tube. Press Sample button, Rinsing will take few minutes. Leaving the distilled water inside the flow cell.
- 4) The screen will display "Turn off the System"
- 5) Turn off the power.

Press key «RINSE» to start or stop the water circulation.

ATTENTION! Cautions

Never leave Sample/Reagent inside the flow cell for a long period of time.

13.7 Troubleshooting table

- | | |
|--|---|
| 4) Non aspiration of sample | <ul style="list-style-type: none"> — Check the pump tube — Check tube pipeline — Rinse and check the flow cell |
| 5) The result of water blank is too high | <ul style="list-style-type: none"> — Check the lamp — Change the distilled water — Rinse the flow cell |
| 6) Repeats are not good | <ul style="list-style-type: none"> — Check if there is air bubble in flow cell |
| 7) Controls out of the target | <ul style="list-style-type: none"> — Reagent validity — Assure the program is correct. — The quality of controls — To test again by changing the mode — Check the flow cell and select a new reagent and control |
| 8) The aspirate volume not stable | <ul style="list-style-type: none"> — Liquid system blocked, Rinse the tube pipeline, or change the tube. |

Appendix: Serial Connection protocol

1. Serial port state:

Data is transmitted at 19200 baud,

With 1 stop bit

8 data bit and Even Parity

2. Message format

- The initial parameter must be the capital letter "R", each line consists of 12 data in following order:

Parameter	Possible values	Space mark
• Initial parameter	"R"	,
• Patient name	Character or NULL	,
• Department	Character or NULL	,
• Sex	Character or NULL	,
	"0" = Male	
	"1" = Female	
• Patient No.	Character	,
• Name of test	Character	,
• Result	Character	,
• Units of measurement	Character	,
• Checksum	Character	

Checksum = All of parameter + All of space mark

3. Send Over

S-2000 sends a character string "TRANSFER FINISH", when the task has been completed.

Example:

R,Biachi,, 0, 43, 1, TP, 136.00, g/L, 0794

R,Rossi,,1, 40, 2, TP, 147.00, g/L, 0765

R,Tony,,0,6,3,TP,35.00,g/L,069d

TRANSFER FINISH

Fig 1.

Remove

Fig 2.

Tighten

14 Warranty

The warranty covers the defective components or faulty assembling of the unit for 1 year since the purchasing date.

Keep purchasing bill as the only valid document for the purchasing date.

The warranty does not cover damages caused for a unit malfunction. The unit must be suitably qualified before being used for valuable samples processing.

Return good are not accepted without previous J.P. SELECTA, s.a.u. approval.

Important note:

Before using the analyzer for the first time:

1. Remove the peristaltic pump protective cover. See Fig 1.
2. Tighten the black screw. See Fig 2.
3. Assembly the cover.

Notice to customers:

The product is made up of various components and various materials that must be recycled or, failing that, deposited in the corresponding debris removal sites when the product's life has been completed or when otherwise it is necessary to dispose of it. To do this, the end user who acquires the product must know the current regulations of each municipality and / or locality based on the waste electrical and electronic equipment. The user who acquires this product must be aware of and responsible for the potential effects of the components on the environment and human health as a result of the presence of hazardous substances. Never place the product in a conventional container of citizen scope if a previous dismantling and knowledge of the components that incorporates. If you do not know the procedure to follow, consult with the city council for more information.